


### Lumitech Incisional Wound Dressing


The Lumitech multi-layer incisional wound dressing is intended to manage post-surgical wounds with low to medium exudate levels. The dressing incorporates **Lumina® Adhesive Technology**, a patented UV switchable adhesive that supports long wear time and a gentle dressing removal. The dressing is compatible with the Prospera Flex® NPWT System.

#### Innovative Features

- Peel and place, sterile multi-layer absorbent design
- Integrated Lumina Adhesive Technology
- Optimized design for application over various incision site locations
- NPWT compatible with dome site positioning flexibility
- Available in a variety of sizes

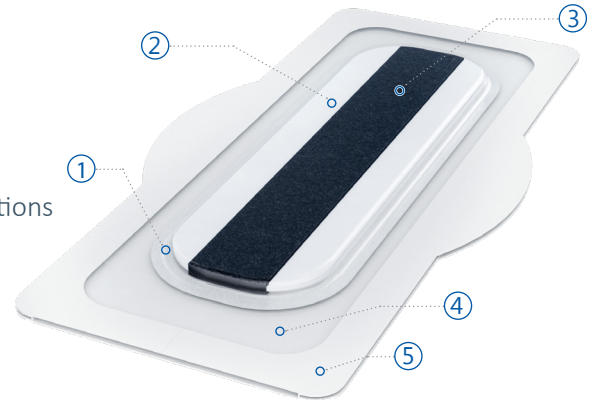

#### LUMITECH®

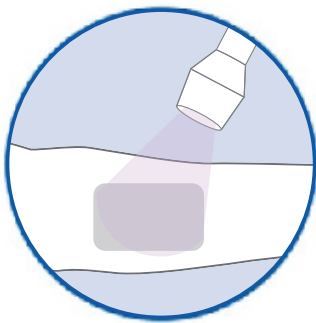


#### Multi Layer Dressing Composition

1. Integrated wound contact layer with extended edges for increased skin integrity and patient comfort
2. Non-woven fibers support up to 20ccs of exudate absorption at the wound site when used independently or with NPWT<sup>1</sup>
3. Polyurethane foam assists with pressure distribution when used with NPWT
4. Polyurethane film incorporating Lumina Adhesive Technology provides flexible, breathable, water-resistant, contaminant blocking protection with a long wear time and gentle removal
5. Release liners and carrier layers for ease of application

#### Lumina Adhesive Technology Benefits<sup>2</sup>

- Adhesive adheres as strong as acrylic dressings, supporting strong seal in unswitched state
- Adhesive releases more gently than silicone following UV light exposure, minimizes pain and periwound trauma at dressing removal

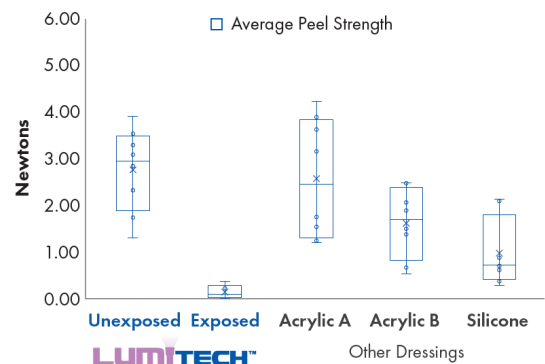

**Compatible with the DeRoyal Prospera Flex NPWT Single-Use System**



Product No.	Description	Packaging
46-LT13	Lumitech 13cm Incisional Dressing	10 Ea/Cs
46-LT20	Lumitech 20cm Incisional Dressing	10 Ea/Cs
46-LTUV	Lumitech Compatible UV Flashlight	10 Ea/Cs
46-LT20S-30	Lumitech Securement Strips (20cm)	6 Ea/Pk, 5 Pks/Cs
46-LT20S-120	Lumitech Securement Strips (20cm)	6 Ea/Pk, 20 Pks/Cs
47-813K-01	Incision Management Convenience Kit (Small)	1 Ea
47-813K	Incision Management Convenience Kit (Small) x1 Prospera Flex NPWT System x2 Lumitech Incisional Wound Dressings (13cm) x1 Lumitech Compatible UV Flashlight Accessories	5 Ea/Cs
47-820K-01	Incision Management Convenience Kit (Medium)	1 Ea
47-820K	Incision Management Convenience Kit (Medium) x1 Prospera Flex NPWT System x2 Lumitech Incisional Wound Dressings (20cm) x1 Lumitech Compatible UV Flashlight Accessories	5 Ea/Cs

**To learn more about our products or to place an order,  
Call 1-888-938-7828 or visit [DeRoyal.com](http://DeRoyal.com).**

1- Data on file, 20cm Dressing

2- Cayce et al. "Light Deactivated Adhesive: The Solution to Maintaining Dressing Integrity and Protecting Skin", Presented at SAWC Spring 2023

DeRoyal, the DeRoyal logo, Prospera Flex, and Lumitech are trademarks or registered trademarks of DeRoyal Industries, Inc. Lumina is a registered trademark of Lumina Adhesives AB.

DeRoyal Industries, Inc. | 200 DeBusk Lane | Powell, TN 37849 | 1.888.938.7828 | DeRoyal.com

© 2025, DeRoyal. All Rights Reserved. | Rev. 2/2025 | LIT0-2608